

Ambu® Perfit ACE® & Mini Perfit ACE®

28 sizes in two collars

DATASHEET

The Ambu Perfit ACE, is a unique adjustable collar. Combined with the Mini Perfit ACE, these two collars offer an incredible 28 sizes ensuring that patients can now have the most effective and accurate fit ever.

The collar is packaged completely flat and has a unique chin-flip feature to ease airway management. Ambu Perfit ACE is a disposable adult adjustable collar for immobilizing the cervical column.

Ambu Perfit ACE is an adjustable collar designed to fit any adult with 16 different settings ranging from neckless to tall sizes.

The Ambu Perfit ACE incorporates ratchet latches that allow the rescuer to size the collar exactly to the patient's neck size. Two safety lock buttons hold the collar in place when adjusted.

Ambu Mini Perfit ACE is an adjustable collar designed to fit small pediatric, pediatric, and small adults with 12 different settings to provide an accurate fit.

* Distortion or shadows may appear on the image depending on the settings of the imaging equipment



KEY BENEFITS

Stores completely flat

CT, MR and X-ray compatible*

Standard size markings on both Adult and Mini

Chin flip feature to ease airway management

Tube hooks are located on the collar body to secure the oxygen tubes

Facilitates intubation

Disposable

When shipped, the default position of the the Adult Ace collar is in the Neckless position, and the default position of the Mini Ace is the Small Pediatric position

Ambu

SPECIFICATIONS

Ambu Perfit ACE

Adjustable sizes	Neckless, Short, Regular and Tall
Number of settings	16 adjustments
Dimensions	554 x 215 x 14 mm (21.8" x 8.5" x 0.6")
Weight	200 g (7 ounces)
Color	White/light grey

Ambu Mini Perfit ACE

Adjustable sizes	Small Pediatric, Pediatric, and Small Adult
Number of settings	12 adjustments
Dimensions	447 x 188 x 14 mm (17.6" x 7.4" x 0.6")
Weight	100 g (3.5 ounces)
Color	Yellow/light grey

MATERIAL


Product	Material
Ambu Perfit ACE	
White/light grey plastic	Polypropylene
Foam	Polyethylene
Ambu Mini Perfit ACE	
Yellow plastic	Polypropylene
Foam	Polyethylene

ORDERING SPECIFICATIONS


Item No.	Description
002 810 001	Perfit ACE Extrication Collar, 20/Case
000 281 107	Mini Perfit ACE Extrication Collar, 30/Case



Ambu® Perfit ACE®

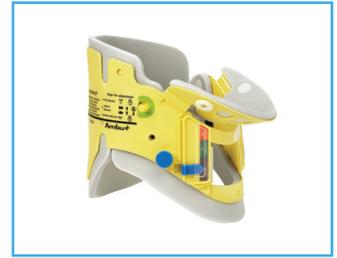

Ambu® Mini Perfit ACE®

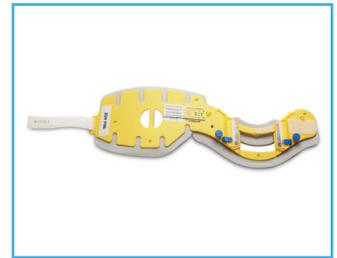

Ambu® Mini Perfit ACE® flat



Ambu USA
6721 Columbia Gateway Drive
Suite 200
Columbia, MD 21045
Tel. 800 262 8462
Fax 800 262 8673
ambuUSA.com

CE US Rx only

Questions? Contact your Ambu representative
at **800-262-8462** or visit **ambuUSA.com**

LBL - 008759 - V01 - 11/2022 - Ambu USA. Technical data may be modified without further notice.