

# FUJI SYSTEMS SOLUTIONS FOR ONE LUNG VENTILATION



**Ambu**

# UNIBLOCKER™


The UNIBLOCKER incorporates a steel mesh polyurethane covered bronchial blocker shaft along with industry leading designs to simplify bronchial blocker placement.

## KEY BENEFITS

- Simple placement for total or segmental bronchial blockade
- Removal of the blocker without detaching the circuit for longterm ventilation
- Placement through standard endotracheal tube

UNIBLOCKER is the only bronchial blocker that incorporates a steel mesh polyurethane covered shaft for direct control of the blocker into the targeted bronchus. This control allows for simple and intuitive placement of the blocker.

In addition to simple placement, the UNIBLOCKER incorporates a gas barrier bronchial cuff which reduces the diffusion of gas during anesthesia.

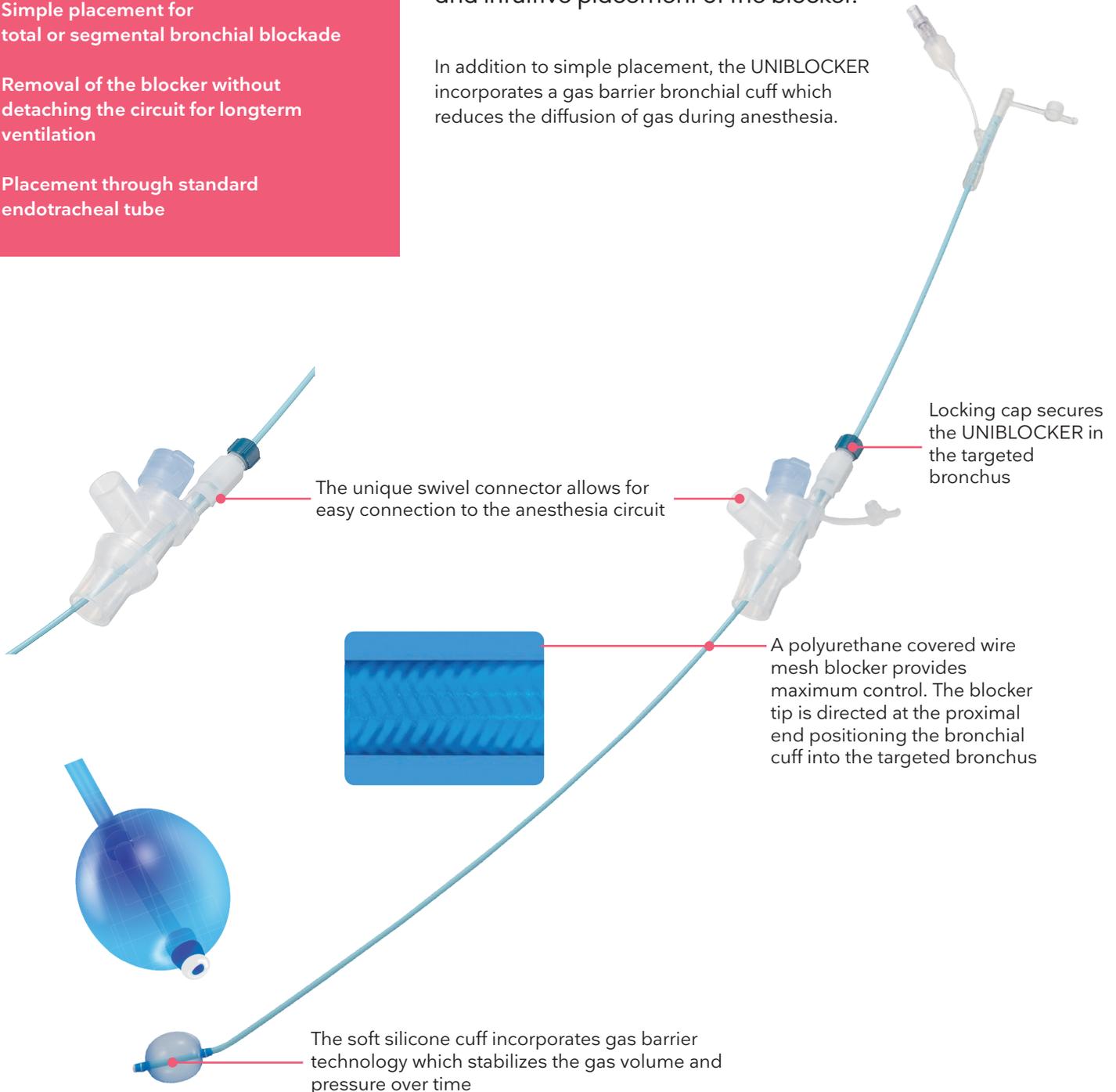

# UNIVENT™

The UNIVENT simplifies bronchial placement by integrating a torque control blocker directly in a single lumen endotracheal tube.


## KEY BENEFITS

- Standard intubation of the UNIVENT into the trachea
- Easy control of the bronchial blocker for left or right bronchial placement
- Polyurethane covered wire mesh bronchial blocker

UNIVENT is a single-lumen silicone endotracheal tube with an integral lumen containing a bronchial blocker.



Simplified placement of the bronchial blocker is achieved with the unique torque control shaft of the bronchial blocker

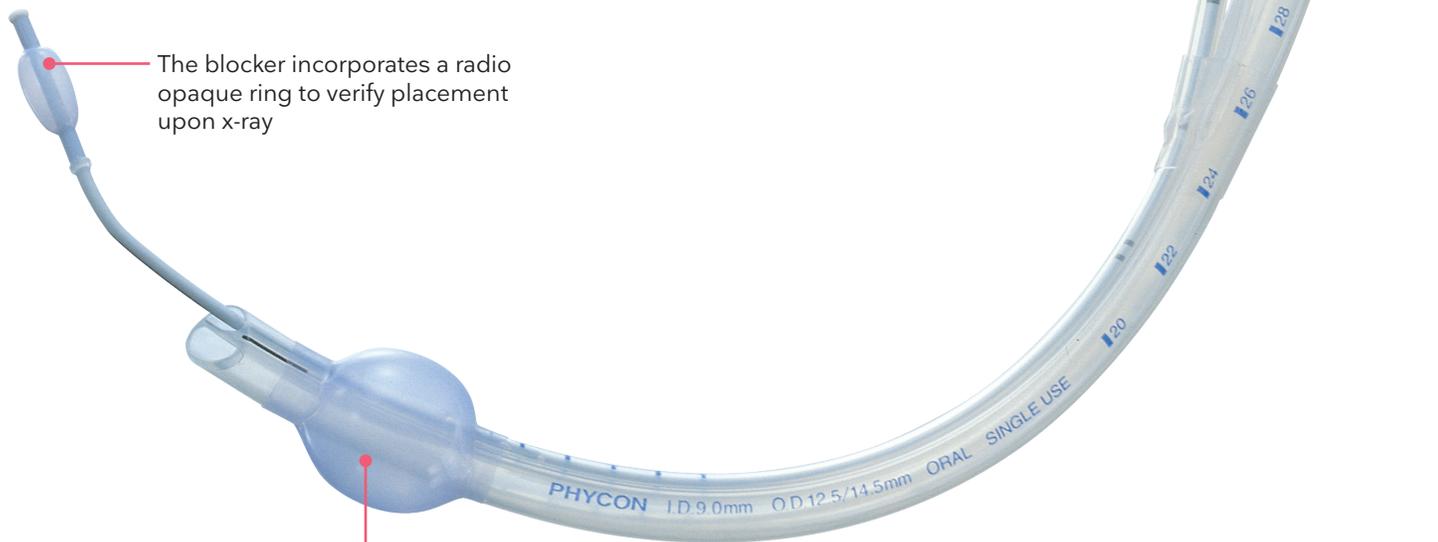

The blocker incorporates a radio opaque ring to verify placement upon x-ray

The blue tracheal high volume low pressure silicone cuff reduces the risk of tear during intubation

# SILBRONCHO™ DOUBLE LUMEN TUBE

The SILBRONCHO double lumen endobronchial tube is designed with medical grade silicone and incorporates a wire reinforced distal and bevel tip to optimize patient safety.

## KEY BENEFITS

- Wire reinforced bronchial tip allows SILBRONCHO to follow the bronchial airway without kinking or dislodging
- The soft flexible tip is designed to reduce the risk of trauma during surgery
- The SILBRONCHO tracheal and bronchial cuffs are less resistant to tear or puncture during intubation

SILBRONCHO utilizes medical grade silicone designed to minimize trauma and reduce heat-related kinking. The wire reinforced distal tip adds flexibility during difficult left or right bronchial placement.



The medical grade silicone is not as thermo sensitive as PVC and reduces the risk of kinking

Wire reinforcement reduces the risk of kinking and provides x-ray verification

Cuff design and location at the distal tip allows for increased depth of placement increasing the margin of safety

The bevel tip reduces the risk of trauma during intubation

The silicone cuff minimizes risk of tearing or puncture during intubation



# SILBRONCHO™ RIGHT SIDED DOUBLE LUMEN TUBE

## KEY BENEFITS

### RIGHT SIDED SILBRONCHO

- Facilitates access to right upper lobe as well as the bronchi of the middle and lower lobes
- The endobronchial cuff is designed so that the right upper lobe can be ventilated without loss of effective lung separation




# WRETT (WIRE REINFORCED ENDOTRACHEAL TUBE)

The Wire Reinforced Endotracheal Tube features a 40 cm, wire reinforced tube to provide adequate depth of placement and ability to seal either the trachea or mainstem bronchus.

## KEY BENEFITS

- Silicone rubber tube provides maximal patient comfort
- Flexible wire-reinforced tube for stability under manipulation and warmth
- Short cuff and longer tube length provide ability to seal in either the trachea or mainstem bronchus

The Wire Reinforced Endotracheal Tube is a single-lumen silicone, wire reinforced tube with a short cuff at the distal end.

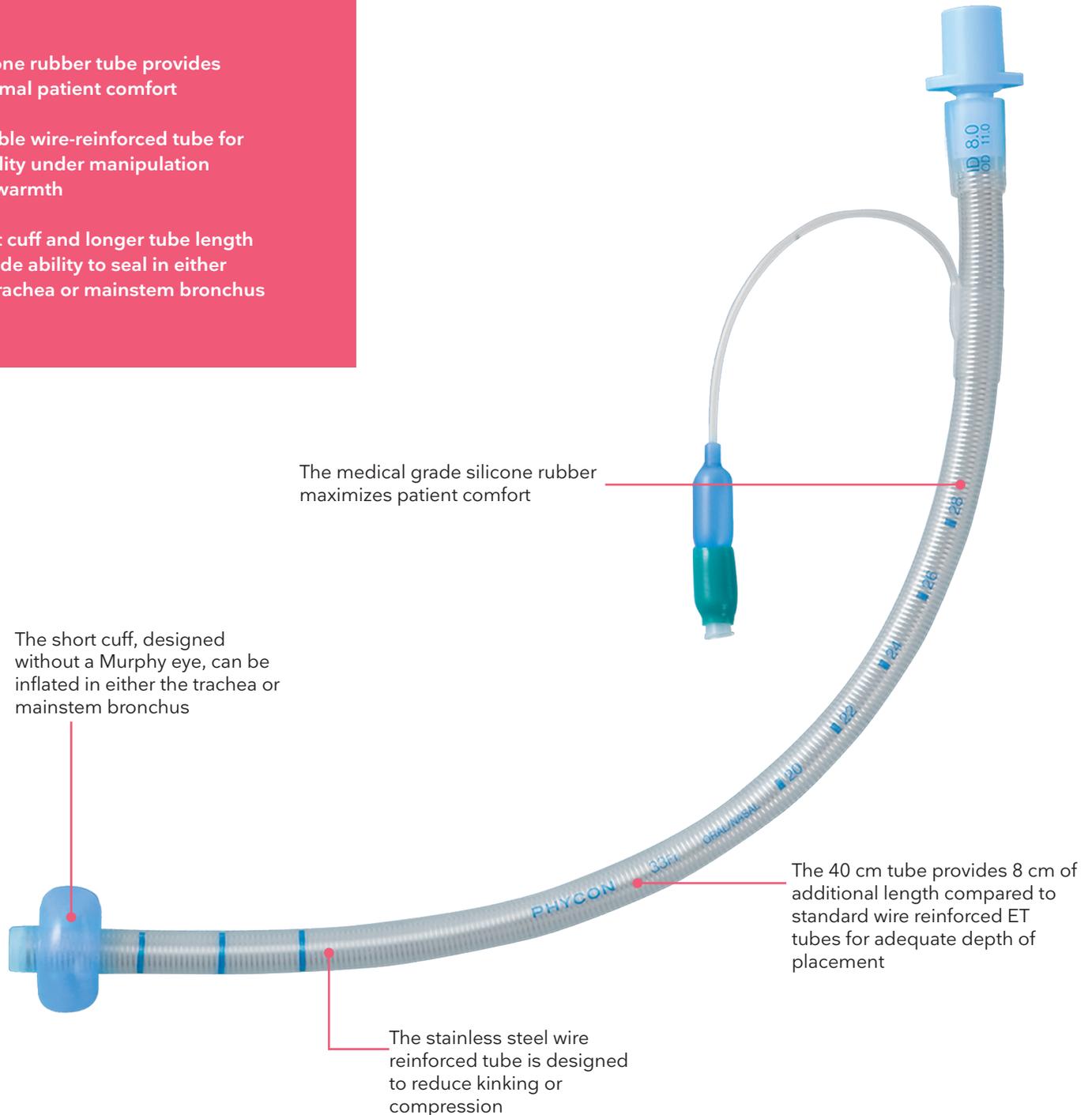

# UNIBLOCKER™

## ORDERING SPECIFICATIONS

Item No.	Description	O.D. (mm)	Total Length (mm)	Maximum Cuff Volume (mm)
1203104	UNIBLOCKER, small	1.7 (5.0 Fr)	505	3
1203105	UNIBLOCKER, adult	3.0 (9.0 Fr)	675	8

# UNIVENT™

## ORDERING SPECIFICATIONS

Item No.	Description	I.D. (mm)	O.D. (mm) (Approx. Value)	Maximum inflation value (Main Cuff) (ml)	Maximum inflation value (Blocker Cuff) (ml)
1202716	UNVENT pediatric without cuff	3.5	7.5/8.0	N/A	3
1202720	UNIVENT pediatric with cuff	4.5	8.5/9.0	15	3
1202927	UNIVENT adult sizes	6.0	9.7/11.5	40	6
1202928	UNIVENT adult sizes	6.5	10.2/12.0	40	6
1202930	UNIVENT adult sizes	7.0	10.7/12.5	40	6
1202931	UNIVENT adult sizes	7.5	11.2/13.0	40	6
1202933	UNIVENT adult sizes	8.0	11.7/13.5	50	6
1202935	UNIVENT adult sizes	8.5	12.2/14.0	50	6
1202936	UNIVENT adult sizes	9.0	12.7/14.5	50	6

# SILBRONCHO™

## ORDERING SPECIFICATIONS

Item No.	Size	I.D. (mm)	O.D. (mm)
1203533	33 Fr left	4.1/6.5	9.5/12.3
1203535	35 Fr left	4.5/7.0	10.0/13.3
1203537	37 Fr left	4.9/7.5	10.5/14.3
1203539	39 Fr left	5.3/8.0	11.0/15.3
1203633	33 Fr right	4.1/6.5	9.5/12.3
1203635	35 Fr right	4.5/7.0	10.0/13.3
1203637	37 Fr right	4.9/7.5	10.5/14.3
1203639	39 Fr right	5.3/8.0	11.0/15.3

# WRETT

## ORDERING SPECIFICATIONS

Item No.	Size	I.D. (mm)	O.D. (mm)
1200824	24 Fr	5.5	8.0
1200828	28 Fr	6.5	9.3
1200831	31 Fr	7.5	10.7

# FUJI SYSTEMS SOLUTIONS FOR ONE LUNG VENTILATION

## Reasons to use the Fuji Systems portfolio for One Lung Ventilation



Flexible high-torque blocker shaft



Unique swivel bronch adapter



Soft, high volume silicone cuff

### Return Goods Policy

The return of goods will only be accepted with the prior agreement in writing by Ambu. When such return of goods is accepted, it is on the understanding that a cancellation or amendment charge may be levied and at no time will the delivery charge be refunded. Such goods must be unopened and in perfect re-saleable condition.

**Ambu** FOREVER  
FORWARD

Ambu, Inc.  
6721 Columbia Gateway Drive, Suite 200  
Columbia, Maryland 21046  
Tel. 800 262 8462  
Fax 800 262 8673  
ambuUSA.com