

DATASHEET

AmbuMan® Basic is designed to fulfill all requirements for the BLS algorithm training.


With the new iQF (Integrated Quality Feedback) App, you can provide training with a high level of precision for CPR performance in adults with cardiac arrest.

CPR's Importance When a Person Suffers from Cardiac Arrest

Brain damage can result if no action is taken within four minutes without oxygen after cardiac arrest. Therefore, it is important that students are trained in this discipline.

iQF, Training That Works for Life

The iQF application includes data-driven information and accurate measurements of relevant CPR parameters. Help your students perform perfect CPR while enjoying games and competing with each other.

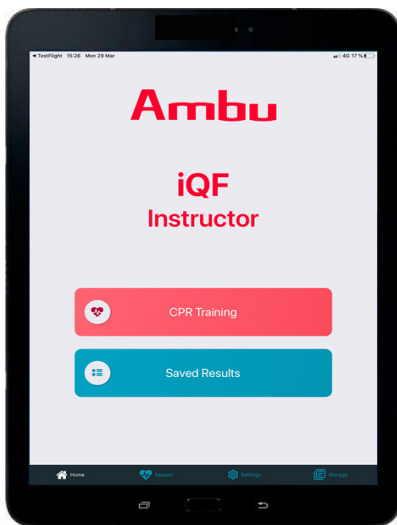

Ambu iQF Instructor

Realistic Ventilation and Correct Compressions

Students will feel how a realistic airway adds resistance and view the rise and fall of the chest when the head is tilted correctly and air is blown into the manikin. The AmbuMan Basic lets students look, listen, and feel the ventilation just like in real life. The unique hygienic system prevents cross-contamination. The correct compression depth, rate, and required ventilation volume can be monitored by the feedback indicator on the manikin.

AmbuMan® Basic iQF



KEY BENEFITS

Robust, reliable, and cost-efficient BLS training manikin with anatomically correct closed torso

Visual feedback built into the manikin torso shows compression depth

Ambu hygienic system for mouth-to-mouth and mouth-to-nose ventilation

Realistic rise and fall of the chest during ventilation

Provides feedback on compression rate and depth and is AHA compliant

TRAIN IQF, TRAIN FOR LIFE

Visualize in detail the compression depth, rate, and ventilation volume at your fingertips with the Instructor and Trainee iQF Apps. Connected by Bluetooth®.

Ambu

SPECIFICATIONS

Weight torso	6 kg
Length	70 cm
Width	37 cm
Weight 4 x torso and bag	30 kg
Max. ventilation volume	1200 ml
Max. compression depth	7 cm

MATERIALS

Basic unit	ABS (Acrylonitrile Butadiene Styrene)
Head	Polyvinyl chloride, hard
Chest skin	Polyvinyl chloride, soft
Face piece	Polyvinyl chloride, soft
Head bag disposable	Polyethylene
Training suit	Polyester
Carrying bag	PVC coated nylon

ORDERING SPECIFICATIONS

Item No.	
A268008000	AmbuMan Basic iQF incl. 1 AmbuMan Basic Torso with training jacket, carry case, 1 face pieces and 25 head bags
234 000 702	Ambu head bags, disposable (set of 100)
234 000 703	Ambu face piece (Set of 5)
234 000 725	Ambu Guedel face piece

* Bluetooth is a registered trademark of Bluetooth SIG, inc.


Questions? Contact your Ambu representative at **800-262-8462** or visit www.ambuusa.com/emergency-care-and-training/training-manikins



Compression Indicator



Unique Ambu Hygienic System

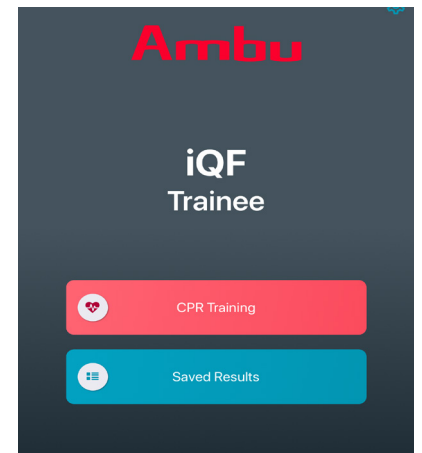

Ambu iQF Trainee



Ambu Inc.
6230 Old Dobbin Lane
Suite 250
Columbia, MD 21045
Tel. 800 262 8462
Fax 800 262 8673
www.ambuUSA.com