

BE MORE PRODUCTIVE WHEREVER YOUR WORK TAKES YOU

Ambu® aScope™ 4 Rhinolaryngo



Ambu



- Perform procedures confidently
- Easily document your procedures
- Streamline your workflow

BE MORE PRODUCTIVE

Whether you are in the outpatient clinic or called to a consult, the Ambu aScope 4 RhinoLaryngo can help you be more productive. This single-use flexible endoscope provides high-quality imaging and connectivity options, making it the perfect grab-and-go solution.

Excellent imaging helps you perform procedures confidently

The aScope 4 RhinoLaryngo is part of a platform that includes the portable, full HD Ambu aView 2 Advance monitor.

You can easily navigate the upper airway and confidently identify anatomical structures with the help of the scope's high-bending angles and precise tip motion combined with superior imaging on the HD monitor.

- Full HD resolution (1920 x 1080) on a 12.8-inch touchscreen provides large, sharp images
- Connect to a high-quality external monitor if a larger image is needed
- Advanced image processing with adaptive light control ensures continuous optimization of the image quality
- Ability to adjust color, contrast, sharpness and brightness settings according to your preferences in varying clinical situations
- Anti-reflective screen surface provides a clear view by reducing ambient light disturbances
- Live 180-degree image rotation to match how you prefer to hold and drive the endoscope

EASILY CONNECT TO MEET YOUR DOCUMENTATION REQUIREMENTS

The aScope 4 RhinoLaryngo endoscope and aView 2 Advance monitor offer a wide range of connectivity options to support the documentation of your procedures. The solution enables you to record images and videos of the procedure, which you can then easily transfer to your Electronic Health Records (EHR) system.



3 SIMPLE STEPS TO DOCUMENTATION

1. **Connect to the hospital network**
using Wi-Fi or Ethernet connection
2. **Record images and video of the procedure**
by pressing record on the aView 2 Advance monitor
3. **Send your recorded patient imaging data**
to the PACS server to enable access from your EHR system



SINGLE-USE EFFICIENCY STREAMLINES YOUR WORKFLOW

Whether you're on a consult or in the clinic, time is of the essence. aScope 4 RhinoLaryngo is always available and easy to take with you streamlining your entire workflow, so you can keep the focus where it counts - on your patients.

Just grab it and go

This lightweight and portable platform is always available to take with you on a consult. With the single-use aScope 4 RhinoLaryngo, you no longer have to wait for a clean scope. The aView 2 Advance monitor is easy to maneuver into a comfortable position close to you and the patient. You can also be prepared for any given clinical situation as the monitor fits into small spaces and rotates 180 degrees to connect from either side. See live images within seconds of powering on the monitor and take advantage of more than three hours of battery life.





Feel the difference of a lightweight, ergonomic design

The lightweight design of the aScope 4 RhinoLaryngo with its ergonomic, rippled-surface handle, will always feel comfortable in your hands, no matter how many procedures you do each day.



Share the view

The aView 2 Advance allows for greater patient involvement. The monitor also records picture and video to allow for post-procedure collaboration.

Benefit from single-use efficiency

The sterile aScope 4 RhinoLaryngo makes it simple to adhere to your cleaning guidelines and improves your productivity. Because it is made to use and discard, you no longer have to worry about the cost and hassle of cleaning and maintaining reusable scopes.



A COMPLETE SOLUTION THAT IMPROVES YOUR EVERYDAY WORKFLOW

The aScope 4 RhinoLaryngo includes a portfolio of two single-use flexible video rhinolaryngoscopes and the high-quality, full HD portable monitor aView 2 Advance:

Ambu® aScope™ 4 RhinoLaryngo Slim

The aScope 4 RhinoLaryngo Slim is an excellent choice for nasal endoscopy and laryngoscopy. It delivers high-quality imaging on the aView 2 Advance monitor, and the small outer diameter helps minimize patient discomfort during procedures. The precise tip motion and high-bending angles allow you to easily maneuver the endoscope in the upper airway.

- Outer diameter: 3.0 mm
- Working channel: n/a
- Length: 300 mm
- Bending angles: 130°/130°



Ambu® aScope™ 4 RhinoLaryngo Intervention

The aScope 4 RhinoLaryngo Intervention is designed for therapeutic procedures. It is ready when you are, requires no post-procedural handling and comes at a low investment cost. Making it particularly suitable for outpatient-clinic biopsies. The Intervention model delivers high-quality imaging and easy maneuverability.

- Outer diameter: 5.0 mm
- Working channel: 2.2 mm
- Length: 350 mm
- Bending angles: 130°/130°



THE SINGLE-USE ENDOSCOPE YOU CAN TRUST

Ambu has the experience, setup, and focus it takes to deliver the single-use solution you need

The world's largest supplier of single-use endoscopes

Since launching the first aScope in 2009, the demand for Ambu aScope technology has steadily risen with increasing numbers of bronchoscopies and intubations being performed worldwide. Bronchoscopes are among the 1.5 million endoscopes that Ambu sold last year. Eighty-seven percent of doctors who evaluated Ambu® aScope™ 4 Broncho in their ICU said it could replace their current reusable scope¹.


Dedicated to single-use endoscopy

We currently have 100 Ambu engineers and 160 strategic partner-engineers working exclusively on single-use endoscopes to make sure you always have the most innovative solutions at your disposal.


Ambu R&D works closely with medical specialists and personnel at our own advanced production facilities to deliver high-quality efficiently. Because we have an agile infrastructure, we are ready to scale up and deliver innovative new products fast.

The Ambu single-use endoscope portfolio:

In addition to our pulmonology and ENT solutions, Ambu has expanded its single-use endoscope portfolio with urology and GI solutions.

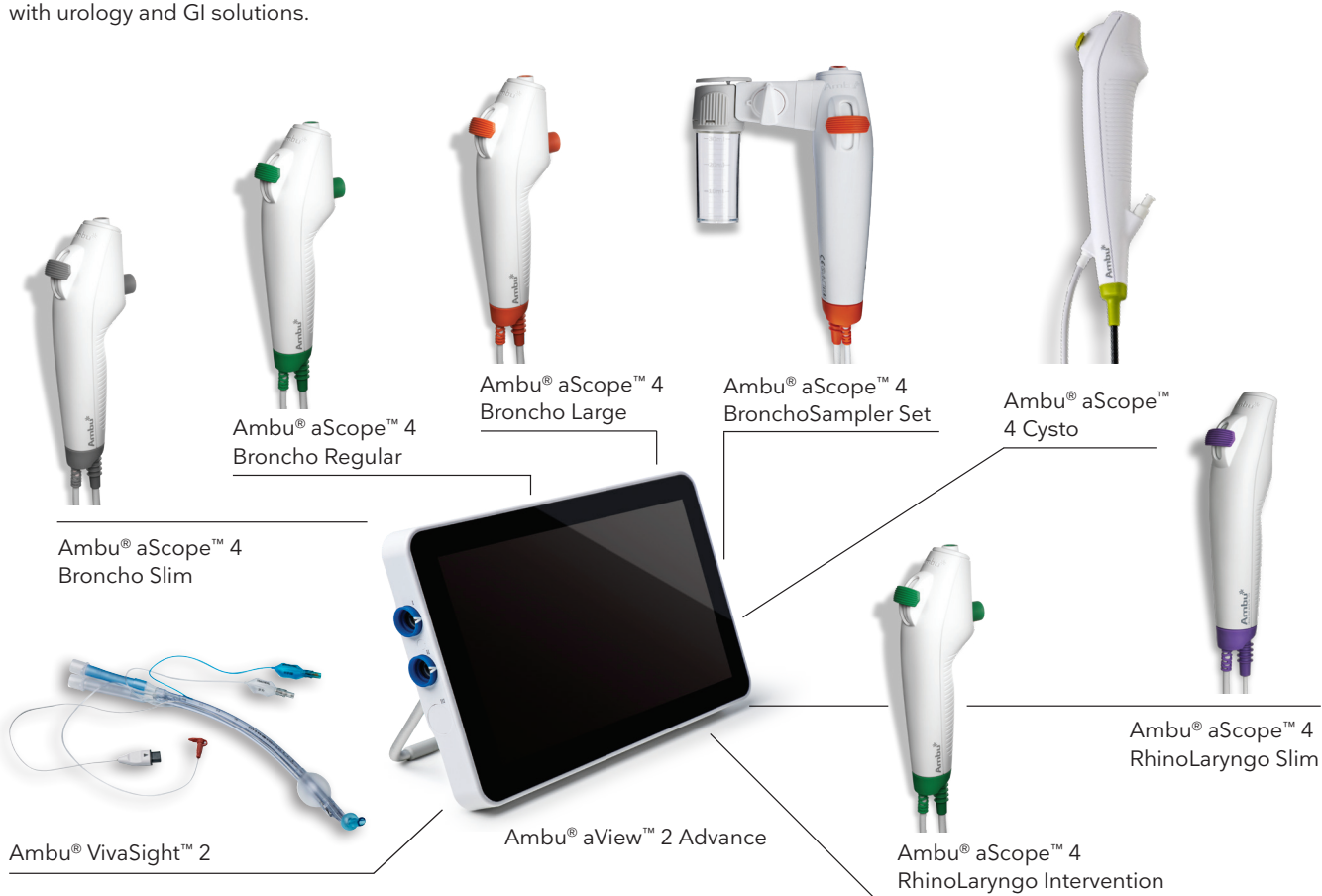

¹JA Livesey and MJG Dunn 2018 "Prospective non-interventional evaluation of intubation and intensive care use of the new Ambu aScope 4 Broncho and aView 2 Advance"

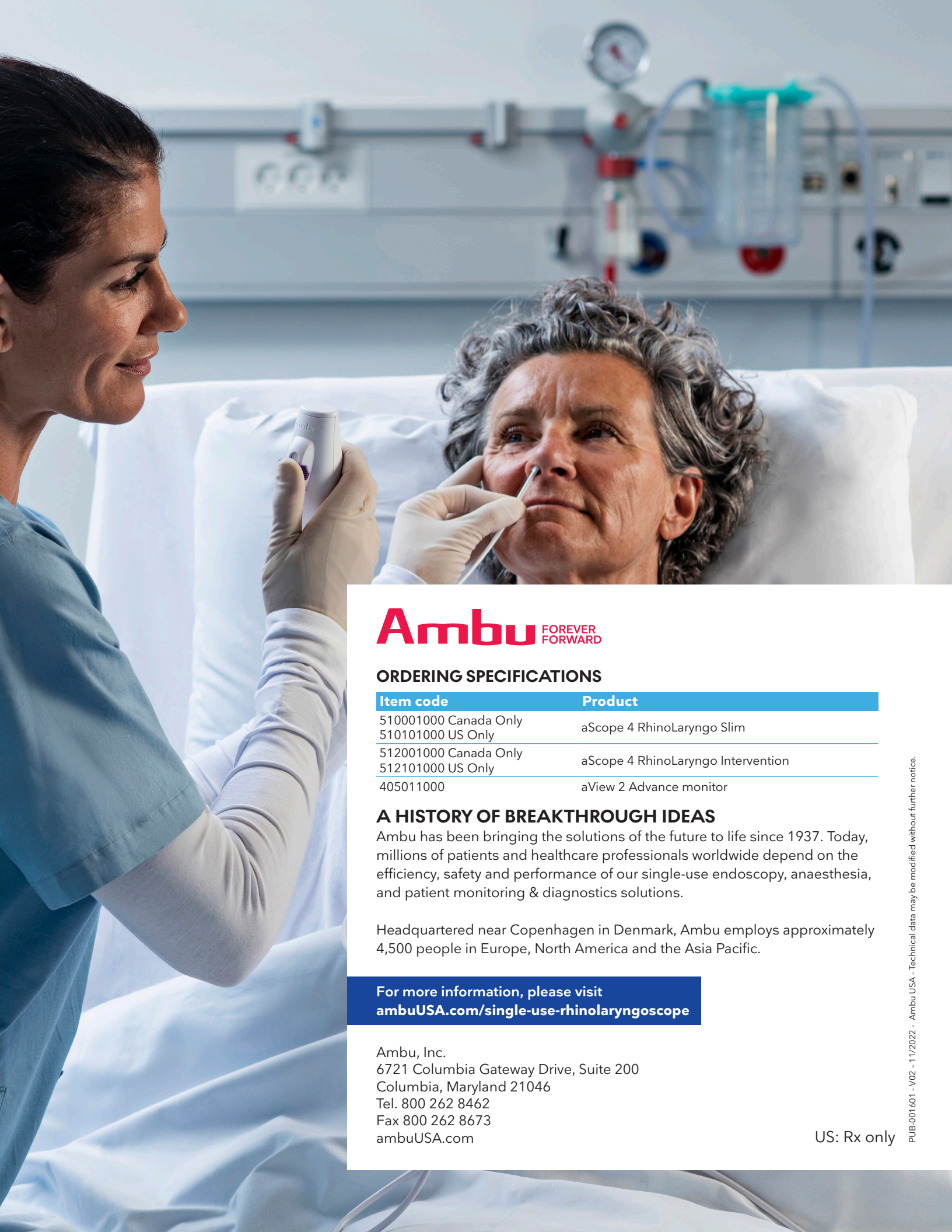

Ambu FOREVER FORWARD

ORDERING SPECIFICATIONS

Item code	Product
510001000 Canada Only 510101000 US Only	aScope 4 RhinoLaryngo Slim
512001000 Canada Only 512101000 US Only	aScope 4 RhinoLaryngo Intervention
405011000	aView 2 Advance monitor

A HISTORY OF BREAKTHROUGH IDEAS

Ambu has been bringing the solutions of the future to life since 1937. Today, millions of patients and healthcare professionals worldwide depend on the efficiency, safety and performance of our single-use endoscopy, anaesthesia, and patient monitoring & diagnostics solutions.

Headquartered near Copenhagen in Denmark, Ambu employs approximately 4,500 people in Europe, North America and the Asia Pacific.

For more information, please visit
ambuUSA.com/single-use-rhinolaryngoscope

Ambu, Inc.
6721 Columbia Gateway Drive, Suite 200
Columbia, Maryland 21046
Tel. 800 262 8462
Fax 800 262 8673
ambuUSA.com

US: Rx only