

Ambu[®] Neuroline 710



Key Benefits

- Special "duck foot" design ensures excellent adhesion
- Flexible soft material
- Repositionable
- Ag/AgCl sensor for best trace quality
- Available in different leadwire configurations

Recommended applications

- Intraoperative Monitoring (IOM)
- Surface Electromyography (EMG)
- Nerve Conduction Studies (NCS)
- Electronystagmography (ENG)
- Electro-Oculography (EOG)
- Polysomnography (PSG)
- Evoked Potentials (EP)
- Reference electrode

Ambu[®] Neuroline 710

Ambu Neuroline 710 is a multiple purpose disposable neurology electrode that can be used for all neurological examinations where the electrodes are placed outside the hair zone. The Ambu Neuroline solid gel electrode is self-adhesive, and it is easy to reposition.

The solid gel electrode is made of a very soft and flexible material, which ensures patient comfort and trace quality even during sleep studies. The 710 electrode is available with white or color-coded lead wires in different lengths.

Specifications

Dimensions	
Electrode size (L x W in mm)	28 x 20
Skin contact size (L x W in mm)	28 x 20
Adhesive area (in mm ²)	490
Height excluding connector/wire (in mm)	1
Sensor	
Sensor material	Silver/silver chloride
Gel system	Solid gel
Sensor area (in mm ²)	490
Gel area/measuring area (in mm ²)	490
Environment	
PVC-free electrode	No
Latex-free electrode	Yes
PVC-free lead wire	No
Latex-free lead wire	Yes
PVC-free packaging	Yes

Materials

Electrode	
Bio-compatible	Yes
Sensor	Silver/silver chloride (Ag/AgCl)
Outer carrier	Carbon filled polyvinyl chloride
Outer adhesive	Solid gel
Upper part	Non-woven polypropylene
Lead-wire core	Flexible copper
Lead-wire insulation	Polyvinyl chloride (PVC)
Lead-wire connector	Brass/polypropylene (PP)
Liner	Siliconized polyester film
Packaging	
Pouches, inner layer	Polyethylene (PE)
Pouches, center layer	Aluminium (Al)
Pouches, outer layer	Paper
Boxes	Cardboard

Ordering Specifications

Item No.	Connector	Lead Wire Length	Qty Per Pack
710 10-K/12	K-K/C	100, 150cm	12
710 10-K/12	K-K/C	150cm	12

Directions

Prepare the skin using an abrasive. Degrease using alcohol swab. Peel the electrode from its protective backing and apply. The electrode can be repositioned several times for better results.

Precautions

Single patient use only. For monitoring and recording biopotential on the skin.

Ambu USA
6230 Old Dobbin Lane
Columbia, MD 21045
Tel. 800 262 8462
Fax 800 262 8673
ambuUSA.com

Connectors




K = 1,5 mm connector



J = 2 mm connector



SC = 0.7 mm connector

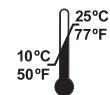

Storage temperature

CE US: Rx only