

Ambu® Silicone Face Mask

DATASHEET

The Ambu® Silicone Face Mask is designed for use with manual and automatic resuscitators, anaesthetic circuits, and ventilators.

Ambu Silicone Face Masks feature a flexible cushion manufactured from 100 % silicone that works to establish a tight yet comfortable seal. The cushion is self-inflating, and the pressure within the cushion may be regulated (except for size 0A).

A transparent dome allows for a visual check of the patient's condition, and an integrated thumb rest supports hand positioning throughout ventilation.

Ambu Silicone Face Masks can be reprocessed up to 30 times. They are suitable for both manual and automated cleaning and sterilization, and are capable of being autoclaved at a temperature of 134 °C (except size 0A). The flexible cushion is also removable to facilitate thorough cleaning.



KEY BENEFITS

Suitable for reprocessing up to 30 times

MR Safe

Transparent dome for visual check of bleeding, vomitus and spontaneous breathing

Thumb rest to support hand positioning for a tight, easy seal

Comfortable design and texture

Flexible cushion manufactured from 100 % silicone which is also removable to allow for comprehensive cleaning

Suitable for autoclaving at a temperature of 134 °C/ 273.2 °F (except size 0A)

Available in 6 sizes

Ambu

MATERIALS


Product	Material
Cuff	Silicone rubber
Cuff size 0A	Silicone rubber with foam filling
Dome	Polysulphone
Mask connector	Silicone rubber

SIZE

Type	Weight
0A	28 g
0 w.boring	39 g
0	38 g
2	63 g
3/4	78 g
4	80 g
5	105 g
6	105 g

ORDERING SPECIFICATIONS

Item No.	Type
000 313 000	Ambu Silicone Face Mask size 0A
000 315 000	Ambu Silicone Face Mask size 0
000 316 000	Ambu Silicone Face Mask size 2
000 312 000	Ambu Silicone Face Mask size 3/4
000 317 000	Ambu Silicone Face Mask size 5
000 319 000	Ambu Silicone Face Mask size 6

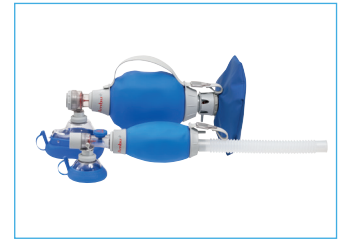

Ambu Silicone Face Mask connected to Ambu Mark IV Resuscitator



Ambu Silicone Face Mask connected to Ambu Oval Plus Silicone Resuscitators

Questions? Contact your Ambu representative at **800-262-8462** or visit **ambuUSA.com**



Ambu, Inc.
 6721 Columbia Gateway Drive
 Suite 200
 Columbia, MD 21046
 Tel. 800 262 8462
 Fax 800 262 8673
 ambuUSA.com

US: Rx only

LBL - 008296 - V01 - 03/2023 - Ambu USA. Technical data may be modified without further notice.